

BOUSSIGNAC

CPAP System

For Acute Pulmonary Edema (APE)

Low O₂ consumption

- Only 15-25 LPM flow achieves CPAP of 5-10 cmH₂O
- "D" cylinder run time 20+ min. at 15 LPM flow

No flow generator

- Boussignac CPAP only needs O₂ source & flowmeter
- No capital equipment to buy and repair
- Completely disposable, single use system
- No down time in the ED waiting to retrieve equipment

Only completely open CPAP system

- Eliminates re-breathing
- Reduces risk of barotrauma/pneumothorax
- Accommodates patient's high peak inspiratory demand
- Permits patient suctioning without removing mask
- Allows use with a nebulizer


Easy to use

- Boussignac CPAP set up in < 2 minutes
- Titrate CPAP simply by adjusting flowmeter
- No valves to change
- Minimal training time

Lightweight and portable

- Complete Boussignac CPAP System weighs only 6.8oz.
- Fits easily into airway bag or pocket (6" x 5" x 4")

Boussignac CPAP works the same way as the turbines of a jet engine.

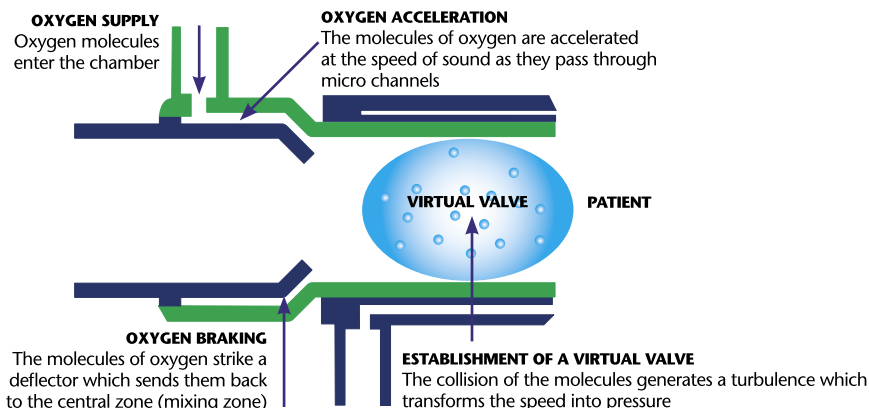

Actual Size

Vitaid
Aiding Life
1.800.267.9301 - www.vitaid.com

BOUSSIGNAC

CPAP System + Nebulizer

For Acute Asthma

30% more drug delivery

- Nebulized medication does not stream to atmosphere
- More effective and cost efficient administration of medicine

Nebulized particle = 1 μ

- Provides for deep lung penetration of nebulized medication

CPAP can be titrated to patient need

- Allows for maintenance of low positive pressure
- Opens bronchioles/alveoli

Reduction of Auto PEEP

- Helps decrease the patient's work of breathing

Open system

- Accommodates patient's high peak inspiratory flow demand
- Reduces risk of barotrauma/pneumothorax



Product code

Item

5562.313	Boussignac CPAP System with size 3 mask (child)
5562.413	Boussignac CPAP System with size 4 mask (small adult)
5562.513	Boussignac CPAP System with size 5 mask (medium adult)
5562.613	Boussignac CPAP System with size 6 mask (large adult)
5570.13	Boussignac CPAP Only – Single use
R-6415V	Nebulizer Kit for Boussignac CPAP – Single use
10-55357	Manometer – Single use
V-51061	15mm/22mm connector for use with ETT/DLT/Trach Tube

Vitaid
Aiding Life

1.800.267.9301 - www.vitaid.com